

Sunflower Warmer

Radiant Warmer



CE 0123



Happier Living Everyday



A BPL Medical Technologies Magazine

Powered by



Scan the product image above with BPL AR App to view the product video of 



Download BPL-AR App to see this catalogue come to life
Download BPL Promise App to experience the latest in medical technology

Simple and Versatile

Total support from postnatal treatment to NICU care.

The most suitable model can be chosen from a cost-effective line up.

Hi-Low/Tilt



Cabinet/Tilt



Stand



Wall Mount



Ideal for
use in
NICU

*Sufficient optional accessories to neonatal care
can be mounted.*




Hi-Low/Tilt

* The drawer shown is
optional

Surgical operations can be performed safely in the NICU

◆ Wide accessible workspace for surgical treatment


- Easy to handle the anesthetic equipment.
- A range of peripheral devices such as resuscitation unit (T-piece resuscitator, oxygen blender and suction unit), phototherapy unit, examination light, weight monitor, syringe pump, etc. can be mounted through F-Rail.

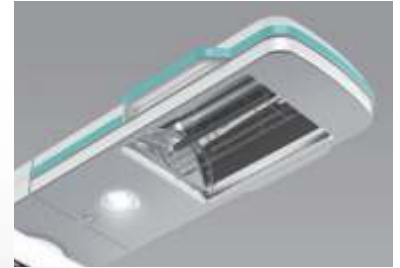

- LED light helps the medical staff to observe the infant. .

“Baby-friendly” function to relieve infant’s stress.


◆ Reduces the infant's stress and supports daily care giving

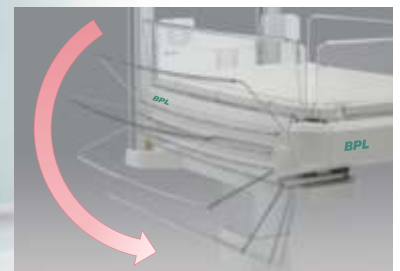

- Baby-guard with very quiet rotary damper.



- Double-wheel caster enables stable and smooth transportation.

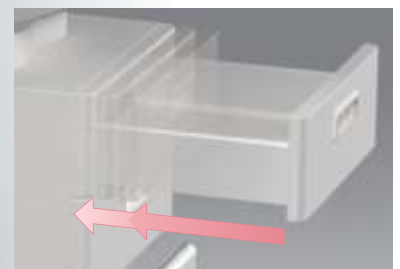

- Self closing doors that close slowly and silently.

*The drawer shown is optional.



Ideal for
Obstetric
use

Best for urgent postnatal treatment.



Hi-Low/Tilt

**The drawer shown is optional*

Prompt and safe resuscitation can be performed.

◆ **Resuscitation treatment in compliance with ILCOR Consensus 2015 can be achieved.**

- CPR Timer prompts the user to assess the health condition of the infant every 30 seconds.



- Resuscitation unit (T-piece resuscitator, oxygen blender and suction unit) and weight monitor, etc. can be mounted to F-Rail.



Utility design.

Sunflower Warmer Hi-Low/Tilt



HL Stand	✓
Fixed Stand	—
Apgar Timer	✓
CPR Timer	✓
Mattress platform tilting	✓
Rotary damper	✓
Drawer	▲
Cabinet	—
Double-wheel Caster	✓
Single-wheel Caster	—
X- Ray cassette tray port	✓
Examination light	✓

✓ : Included as standard
— : Not available ▲ : Optional

Sunflower Warmer Cabinet/Tilt



HL Stand	—
Fixed Stand	✓
Apgar Timer	✓
CPR Timer	✓
Mattress platform tilting	✓
Rotary damper	✓
Drawer	—
Cabinet	✓
Double-wheel Caster	✓
Single-wheel Caster	—
X- Ray cassette tray port	✓
Examination light	✓

✓ : Included as standard — : Not available

Maneuverable design

Sunflower Warmer Stand



HL Stand	—
Fixed Stand	✓
Apgar Timer	✓
CPR Timer	✓
Mattress platform tilting	—
Rotary damper	—
Drawer	—
Cabinet	—
Double-wheel Caster	—
Single-wheel Caster	✓
X- Ray cassette tray port	—
Examination light*	✓

✓ : Included as standard — : Not available

* Except for US

Space-saving design.

Sunflower Warmer Wall Mount



HL Stand	—
Fixed Stand	—
Apgar Timer	✓
CPR Timer	✓
Mattress platform tilting	—
Rotary damper	—
Drawer	—
Cabinet	—
Double-wheel Caster	—
Single-wheel Caster	—
X- Ray cassette tray port	—
Examination light*	✓

✓ : Included as standard — : Not available

* Except for US





Drawer, 2 Layers

For Infa Warmer i/
Sunflower Warmer

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

Self-closing drawer that closes slowly and silently



Handle

For F-Rail

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

The vertical position of the handle can be adjusted easily according to the height of the operator. It is conveniently used to move Sunflower Warmer to another place.



Drainage Hanger

For Infa Warmer i/
Sunflower Warmer

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

The Drainage Hanger is used to hang a low-pressure continuous aspirator and drainage bag on Sunflower Warmer.



Tube Holder

For Infa Warmer i/
Sunflower Warmer/V-505

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

The Tube Holder is attached to the baby guard at the head end of the warmer to hold the patient circuit securely in place.



Shelf

For Sunflower Warmer

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

Shelf with wide space can place a SpO2 monitor, etc.



Tray Set (S)

For F-Rail

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

The vertical position of the Tray Set (S) can be adjusted easily according to the height of the operator. Small instruments necessary for routine medical procedures can be placed on this small tray.



Basket

w/arm For F-Rail

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

The vertical position of the basket can be adjusted easily as desired. The basket is deep enough to conveniently accommodate tall containers such as bottles of alcohol disinfectant solutions.



Monitor Table for Patient monitor

For F-Rail

Hi-Low/Tilt Cabinet/Tilt Stand
Wall Mount

Patient monitor can be mounted on the monitor table



MF-Rail

For F-Rail
(Common to Left & Right)

Hi-Low/Tilt **Cabinet/Tilt** **Stand**
Wall Mount

This rail is used to mount the peripheral devices attachable to the MF-Rail for F-Rail.



I.V Pole

For F-Rail

Hi-Low/Tilt **Cabinet/Tilt** **Stand**
Wall Mount

With the load capacity of 10kg, this I.V Pole can mount the weight of various peripheral devices



Resusci Flow

Hi-Low/Tilt **Cabinet/Tilt** **Stand**
Wall Mount

"Resusci Flow" is a resuscitation unit that performs artificial respiration with a stable pressure.

- Simple operation supports resuscitation.
- Safety function to protect newborn's respiration.
- Supports artificial respiration in various situations.



Infalite

Hi-Low/Tilt **Cabinet/Tilt** **Stand**
Wall Mount

With a LED light and hands-free motion sensor, both brightness and hygiene have been considered with this examination light.

*Either F-Rail attachment or Pole attachment is required for the installation of InfaLite.



BILI-THERAPY Spot Type Arm Type

Hi-Low/Tilt **Cabinet/Tilt** **Stand**
Wall Mount

Effective Phototherapy Treatment and Space-saving

- Effective phototherapy treatment by blue LEDs
Promote early hospital discharge by reducing the treatment time in providing high irradiation intensity and effective wavelength.
- Space-saving
The unit can be mounted on the F-Rail without using a stand for space-saving around the incubator or warmer.
Switchable to examination lamp with one-touch of a button.
The condition of an infant can be monitored by the examination lamp (White LEDs).
The dimmer function ideal for developmental care is available.

Products Lineup

Sunflower Warmer

Type	Description	Type	Description
Hi-Low/Tilt	AC120V±10%	Stand	AC120V±10%
	AC230V±10%		AC230V±10%
Cabinet/Tilt	AC120V±10%	Wall Mount	AC120V±10%
	AC230V±10%		AC230V±10%

Specifications


Power requirements	Customer specified		
Power consumption	Main Unit : 700VA HL Stand : 120VA		
Frequency	50/60Hz		
Type of protection	Class I		
Degree of protection	Type BF applied part		
Classification	IIb, Rule 9		
Heater capacity	500W		
Dimensions	HL/Tilt : 640(W)x1140(D)x1770~1970(H) mm Cabinet/Tilt : 610(W)x1060(D)x1870(H) mm Stand : 730(W)x950(D)x1840(H) mm Wall mount : 450(W)x950(D)x300(H) mm		
Mattress surface	HL/Tilt	800~1000 mm	
	Cabinet/Tilt	900 mm	
	Stand	660~860 mm	
	Wall mount	Depends on mount height	
Weight	HL/Tilt	Approx. 89Kg	
	Cabinet/Tilt	Approx. 96Kg	
	Stand	Approx. 51Kg	
	Wall mount	Approx. 20Kg	

Accessories	HL/Tilt	Skin temp. probe x 1, Dust cover x 1, Cord holder x 2, Operation Manual x 1
	Cabinet/Tilt	Skin temp. probe x 1, Dust cover x 1, Cord holder x 2, Operation Manual x 1
	Stand	Skin temp. probe x 1, Cord holder x 2, Operation Manual x 1
	Wall mount	Skin temp. probe x 1, Operation Manual x 1
Heater output setting range		0~100% (in 5% increments)
Heater output indication		In 10 levels
Skin temp. setting range		34.0~38.0℃ / 93.2~100.4℉ (in 0.1℃/℉ increments)
Skin temp. display range		30.0~42.0℃ / 86.0~107.6℉ (in 0.1℃/℉ increments)
Skin temp. display accuracy		±0.3℃ / ±0.5℉
Timer		0~60 minutes 00 second (in one second increments) Apger mode / CPR mode selectable
Mattress platform tilting		Longitudinal tilt of ±13°(freely adjustable)
Canopy mobility range		Horizontal rotation : 170°(85°both sides)
Lighting		600Lx
Preheat		Available
Alarms		Set temp., Baby check, Power failure, System failure, Skin temp. probe


Dimensional drawing

Unit:mm


■ Hi-Low/Tilt

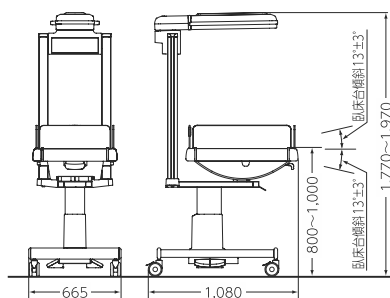



■ Cabinet/Tilt

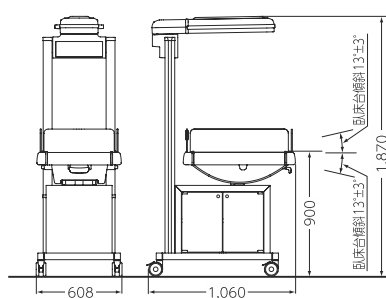

■ Stand

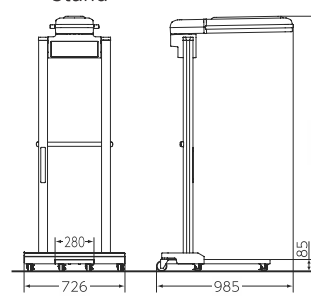

■ Wall Mount

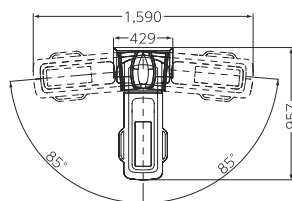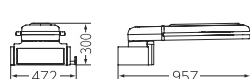

Atom Medical products are produced under strict production control in Japan.
We provide value and trust of the highest quality for the tiny baby lives in the world.

Manufacturer: ATOM MEDICAL CORPORATION

18-15, Hongo 3-chome, Bunkyo-ku, Tokyo, 113-0033 JAPAN
<http://www.atomed.co.jp>

Exporter: ATOM MEDICAL INTERNATIONAL, INC.

Iwakata Bldg., 18-16, Hongo 3-chome, Bunkyo-ku, Tokyo, 113-0033 JAPAN
Phone: + 81 3 3815 2941 Fax: + 81 3 3812 9670/4080

BPL MEDICAL TECHNOLOGIES PRIVATE LIMITED

11th KM, Bannerghatta Road, Arakere, Bangalore - 560076
CIN U33110KA2012PTC067282
Phone +91 80 26484388/ 2648 4350
Toll Free 1800-425-2355
Website www.bplmedicaltechnologies.com
For Enquiries: sales.medical@bpl.in
Copyright©2017 BPL Medical Technologies Private Limited.
All rights reserved.