



## Flat Top Medical Treatment Table



### OPTIONAL



Easy Mobility



Walk free switch

OJAS is a flat top medical treatment table that can be used for general examination and massage treatments. The height can be electronically adjusted through a linear actuator from Linak DENMARK. This table is supplied with an external, detachable & adjustable, face cradle that can be fitted on both sides of the table. The treatment table has a strong powder coated metal base and top mattress is 40 density cushion upholstered with soft to touch premium leatherette.

**TOP:**

Top Design	Flat Top with adjustable face cradle attachment
Top Cushion Measurements	1995 x 750 mm
Foam	50 mm thick 40 Density Foam
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame
Adjustment (Back Rest)	Flat Top Design
Adjustment (Height)	Electronic (Linak Motor)

**BASE & FRAME:**

Base	Powder Coated Metal Frame
Minimum Height	530 mm
Maximum Height	890 mm
Lifting Capacity	Upto 360 Kgs

**UPHOLSTERY:**

Type	Flame/UV/Abrasion Resistant 1.2 mm Coated Fabric
Available Colour Options	25 colour options
Weight	660 gsm + 5%
Abrasion	EN ISO 5470=2, 25000 Cycles

**ACCESSORIES:**



Detachable External Face Cradle
---------------------------------

**GUARANTEES:**

Structure	1 year
Electronic Motors & Controls	3 years


**WEIGHTS:**

Table Weight	84kgs
Shipping Weight	135kgs

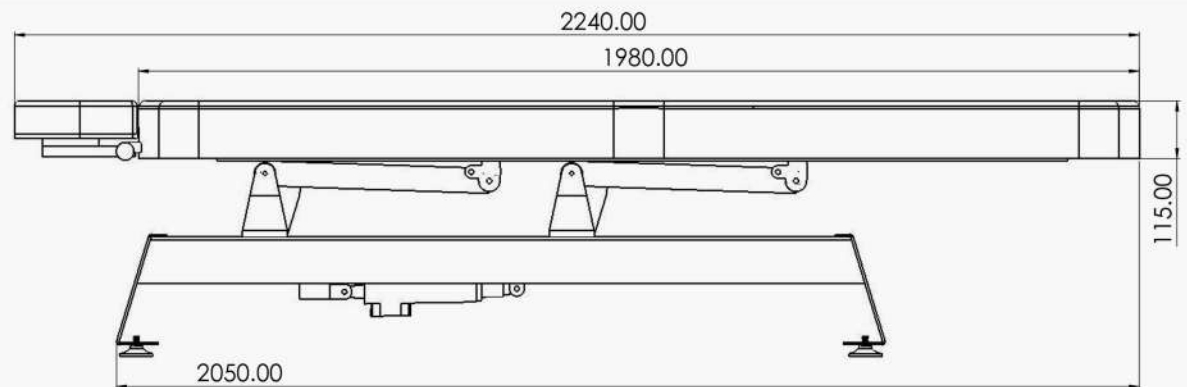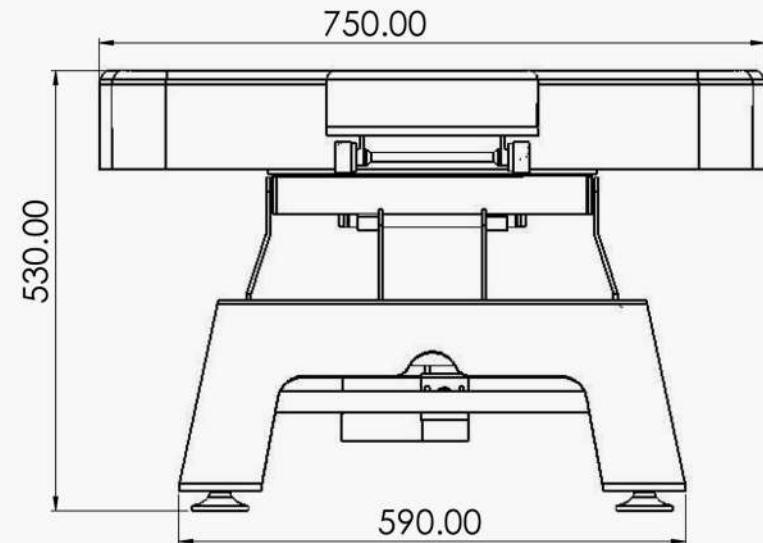


All dimensions are in mm

## MINIMUM HEIGHT

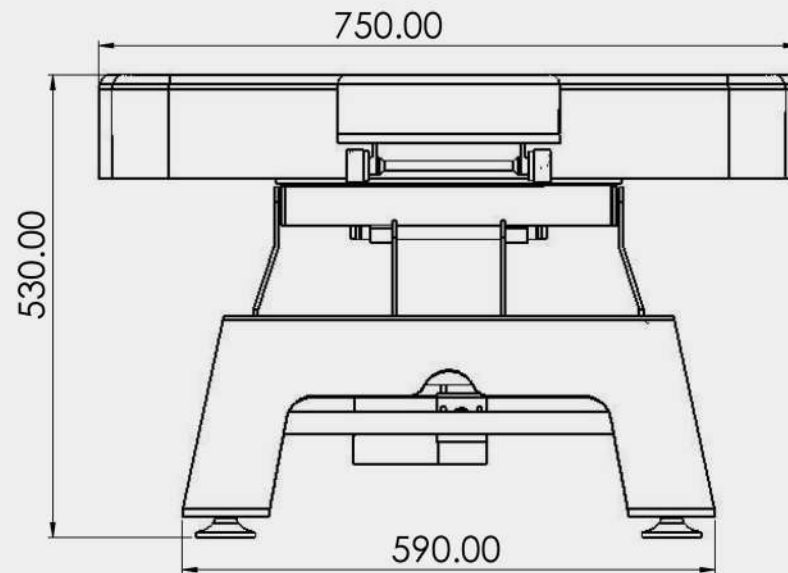

## MAXIMUM HEIGHT

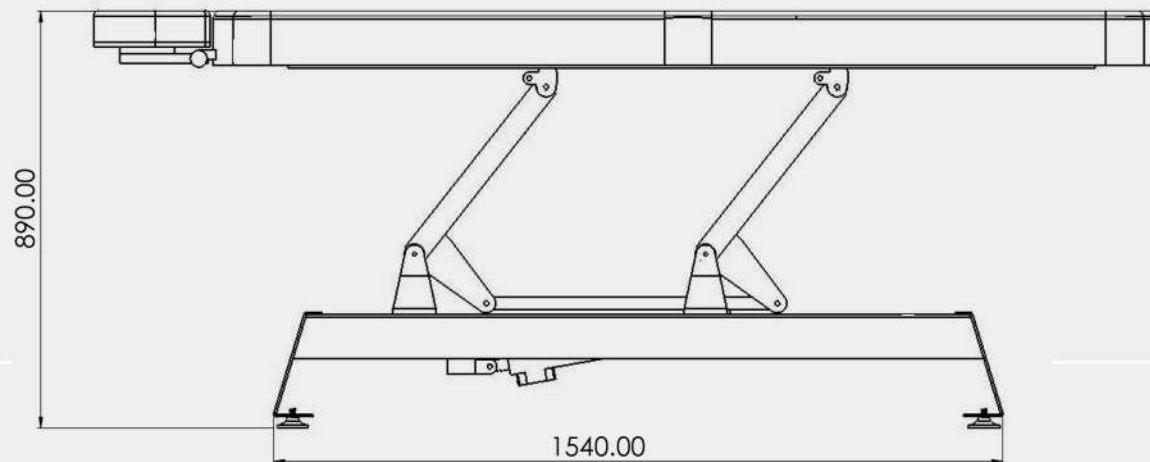

All dimensions are in mm

## 2 Section Medical Treatment Table



Nirav is two section medical treatment table with electronic height and manual headrest adjustment. The headrest position can be adjusted and locked at various positions by a pneumatic gas spring ; whereas the height is controlled through a linear actuator from Linak DENMARK. There is a built in nose slot for breathing with plug cushion. The treatment table has a strong powder coated metal base and top mattress is 40 density cushion upholstered with soft to touch premium leatherette.

### OPTIONAL



Easy Mobility



Walk free switch

# NIRAV





**TOP:**

Top Design	Two Piece articulating top with headrest incline and decline
Top Cushion Measurements	1995 x 750 mm
Foam	50 mm thick 40 Density Foam
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame
Adjustment (Back Rest)	Manual ( Gas Spring )
Adjustment (Height)	Electronic (Linak Motor)

**BASE & FRAME:**

Base	Powder Coated Metal Frame
Minimum Height	620 mm
Maximum Height	890 mm
Lifting Capacity	Upto 360 Kgs

**UPHOLSTERY:**

Type	Flame/UV/Abrasion Resistant 1.2 mm Coated Fabric
Available Colour Options	25 colour options
Weight	660 gsm + 5%
Abrasion	EN ISO 5470=2, 25000 Cycles

**ACCESSORIES:**



Available Colour Options	25 colour options
--------------------------	-------------------

**GUARANTEES:**

Structure	1 year
Electronic Motors & Controls	3 years


**WEIGHTS:**

Table Weight	84kgs
Shipping Weight	135kgs

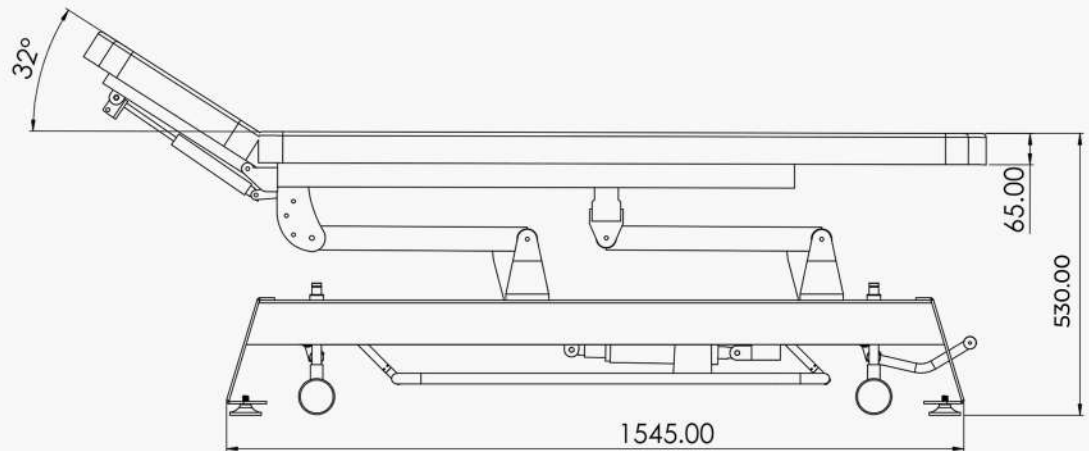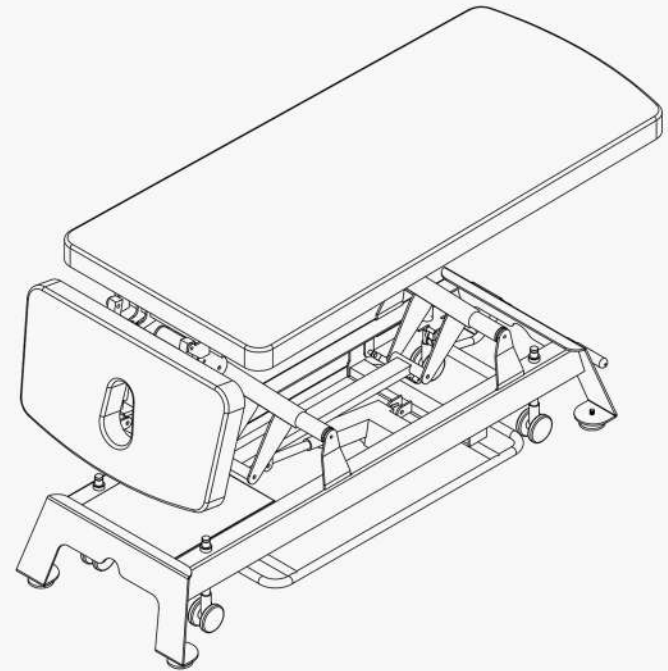


All dimensions are in mm

## MINIMUM HEIGHT

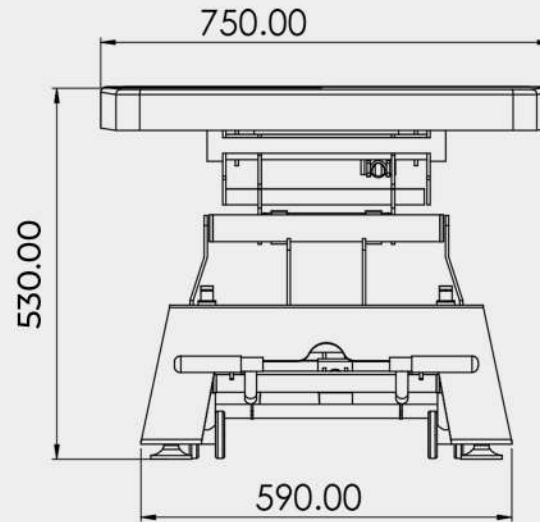

## MAXIMUM HEIGHT

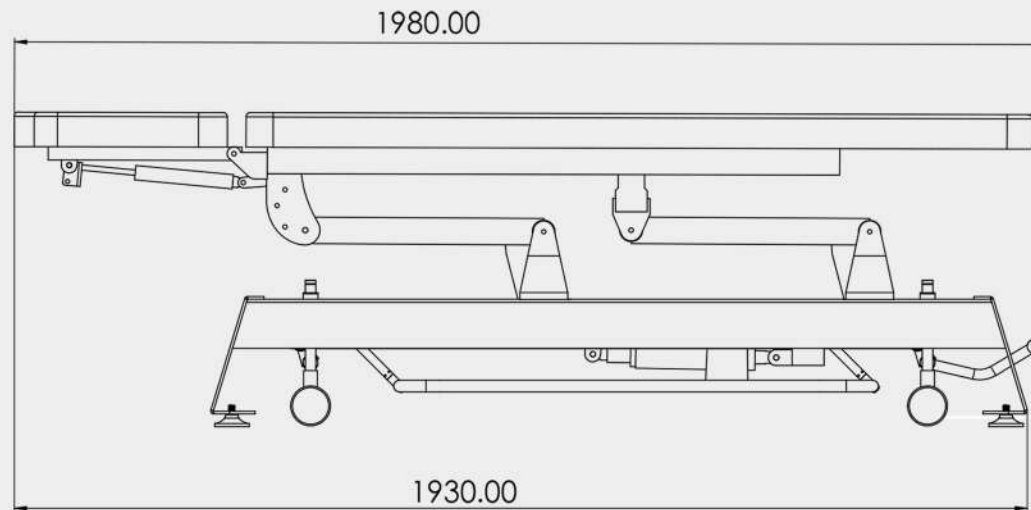


All dimensions are in mm

## 2 Section Medical Treatment cum Examination Table



Manan is our elegant in design and functional two section medical treatment cum examination table with electronic height and manual backrest adjustment. The backrest can be adjusted and locked at various positions by a pneumatic gas spring; whereas the height is controlled through a linear actuator from Linak DENMARK. There is a built-in nose slot for breathing with a plug cushion. The treatment table has a strong powder-coated metal base and the top mattress is 40 density cushion upholstered with soft to touch premium leatherette.

### OPTIONAL

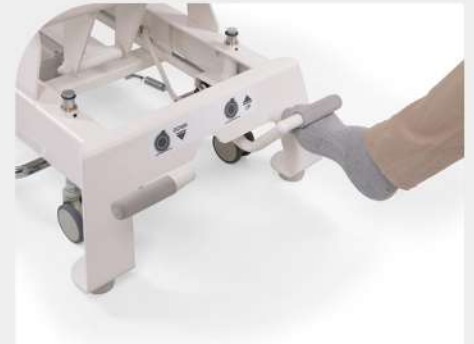

Easy Mobility



Walk free switch



# MANAN



**TOP:**

Top Design	Two Piece articulating top with backrest incline feature
Top Cushion Measurements	2050 x 750 mm
Foam	50 mm thick 40 Density Foam
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame
Adjustment (Back Rest)	Manual ( Gas Spring )
Adjustment (Height)	Electronic (Linak Motor)

**BASE & FRAME:**

Base	Powder Coated Metal Frame
Minimum Height	530 mm
Maximum Height	890 mm
Lifting Capacity	Upto 360 Kgs

**UPHOLSTERY:**

Type	Flame/UV/Abrasion Resistant 1.2 mm Coated Fabric
Available Colour Options	25 colour options
Weight	660 gsm + 5%
Abrasion	EN ISO 5470=2, 25000 Cycles

**ACCESSORIES:**



Available Colour Options	25 colour options
--------------------------	-------------------

**GUARANTEES:**

Structure	1 year
Electronic Motors & Controls	3 years


**WEIGHTS:**

Table Weight	84kgs
Shipping Weight	135kgs

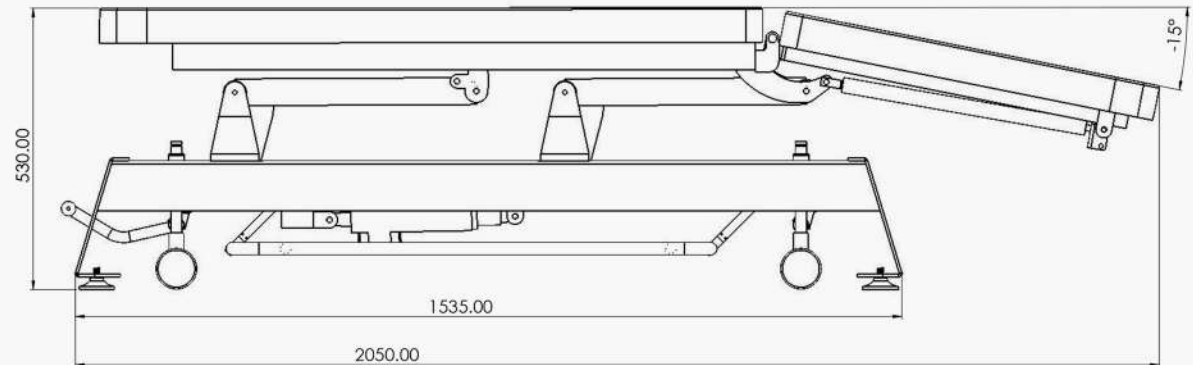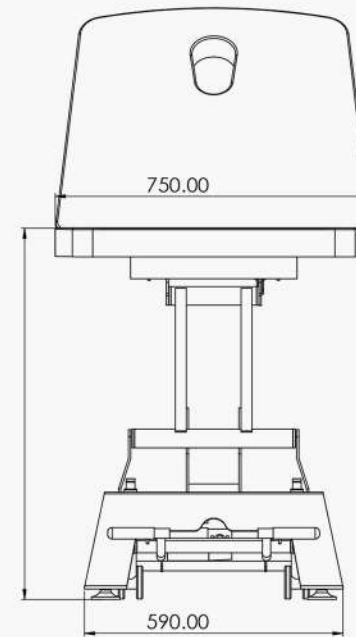


All dimensions are in mm

## MINIMUM HEIGHT

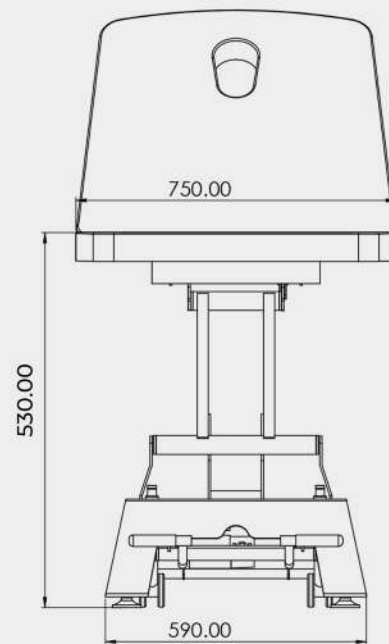

## MAXIMUM HEIGHT

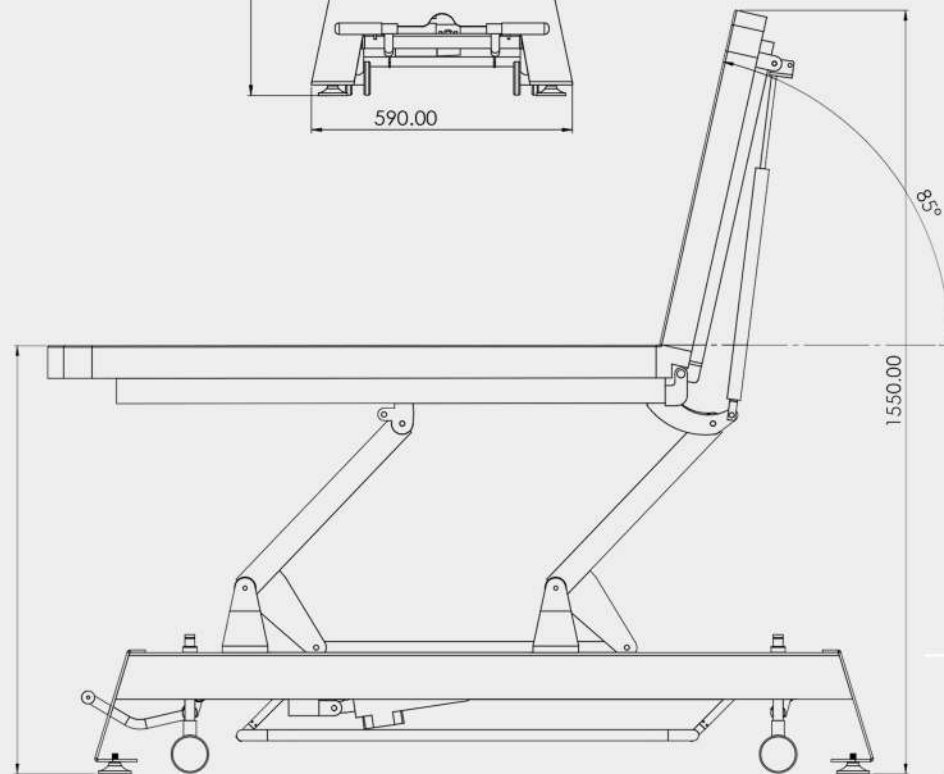

All dimensions are in mm

## 2 Section Physiotherapy Table



Parth is best suited for physiotherapy and physical therapies. Two fold treatment table with electronic height, manual backrest and pneumatically lowerable armrests. The backrest and the armrests can be adjusted and locked at various positions by a pneumatic gas spring; whereas the height is controlled through a linear actuator from Linak DENMARK. Best suited for doctors and physiotherapists to work on the back of the patient. There is a built in nose slot for breathing with plug cushion. The treatment table has a strong powder coated metal base and top mattress is 40 density cushion upholstered with soft to touch premium leatherette.

### OPTIONAL



Easy Mobility



Walk free switch



Adjustable Head Rest is optional

# PARTH

01



02



03



04





**TOP:**

Top Design	Two Piece articulating top with lowerable armrests
Top Cushion Measurements	1995 x 750 mm
Foam	50 mm thick 40 Density Foam
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame
Adjustment (Back Rest)	Manual ( Gas Spring )
Adjustment (Height)	Electronic ( Linak Motor )

**BASE & FRAME:**

Base	Powder Coated Metal Frame
Minimum Height	530 mm
Maximum Height	890 mm
Lifting Capacity	Upto 360 Kgs

**UPHOLSTERY:**

Type	Flame/UV/Abrasion Resistant 1.2 mm Coated Fabric
Available Colour Options	25 colour options
Weight	660 gsm + 5%
Abrasion	EN ISO 5470=2, 25000 Cycles

**ACCESSORIES:**



Available Colour Options	25 colour options
--------------------------	-------------------

**GUARANTEES:**

Structure	1 year
Electronic Motors & Controls	3 years


**WEIGHTS:**

Table Weight	84kgs
Shipping Weight	135kgs

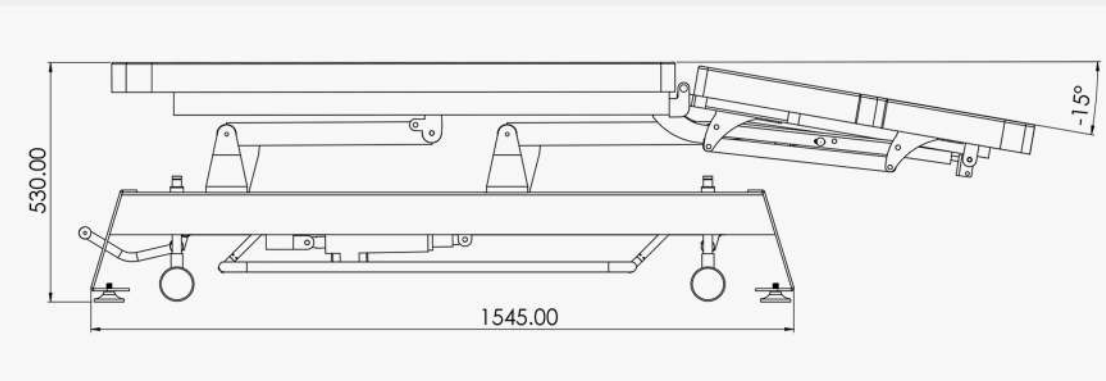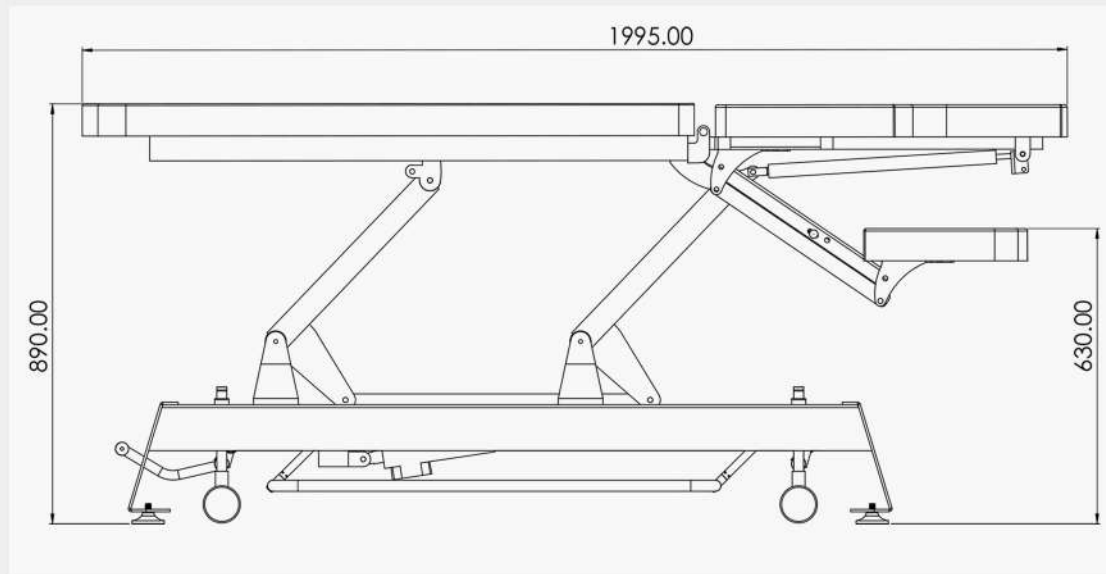


All dimensions are in mm

## MINIMUM HEIGHT

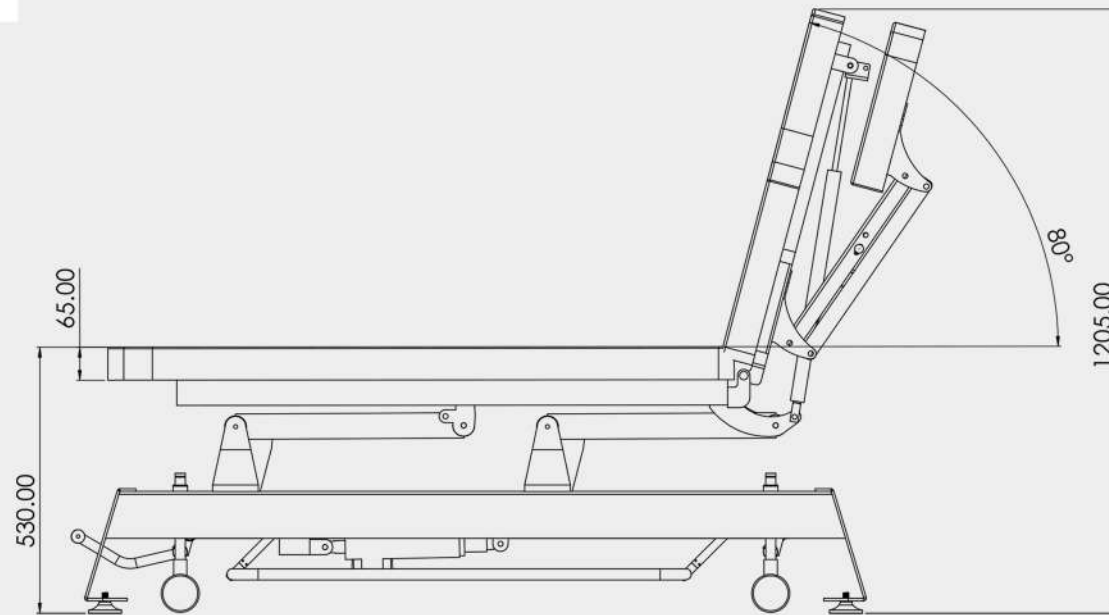

## MAXIMUM HEIGHT

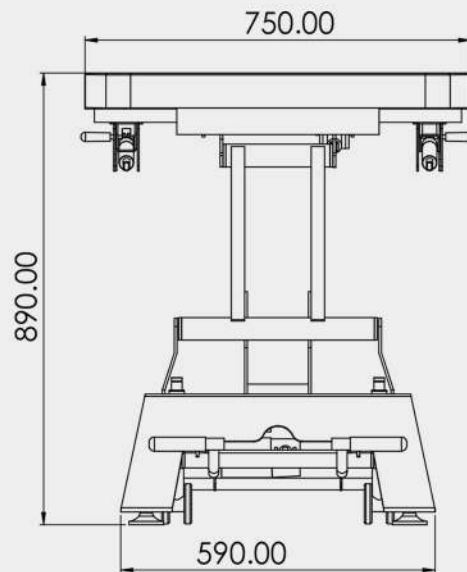

All dimensions are in mm

# MOKSHIT

## 3 Fold Treatment Table



MOKSHIT is a 3 fold treatment table with a variety of functions. It can be used as a general examination table with independently adjustable backrest and leg rest section. Both leg and back sections are controlled through a quality gas spring. You can make it a flat top bed or bring it into a chair position. The height can be electronically adjusted through a linear actuator from Linak DENMARK. There is built in nose slot for easy breathing and plug cushion. Strong powder coated metal base and top mattress is 40 density cushion upholstered with soft to touch premium leatherette.

### OPTIONAL



Easy Mobility



Walk free switch

# MOKSHIT



**TOP:**

Top Design	Three Piece articulating top with adjustable backrest and leg rest position
Top Cushion Measurements	2010 x 750 mm
Foam	50 mm thick 40 Density Foam
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame
Adjustment (Back Rest)	Manual ( Gas Spring )
Adjustment (Height)	Electronic (Linak Motor)
Adjustment (Leg Rest)	Manual

**BASE & FRAME:**

Base	Powder Coated Metal Frame
Minimum Height	525 mm
Maximum Height	890 mm
Lifting Capacity	Upto 360 Kgs

**UPHOLSTERY:**

Type	Flame/UV/Abrasion Resistant 1.2 mm Coated Fabric
Available Colour Options	25 colour options
Weight	660 gsm + 5%
Abrasion	EN ISO 5470=2, 25000 Cycles

**ACCESSORIES:**



Available Colour Options	25 colour options
--------------------------	-------------------

**GUARANTEES:**

Structure	1 year
Electronic Motors & Controls	3 years

**WEIGHTS:**


Table Weight	84kgs
Shipping Weight	135kgs

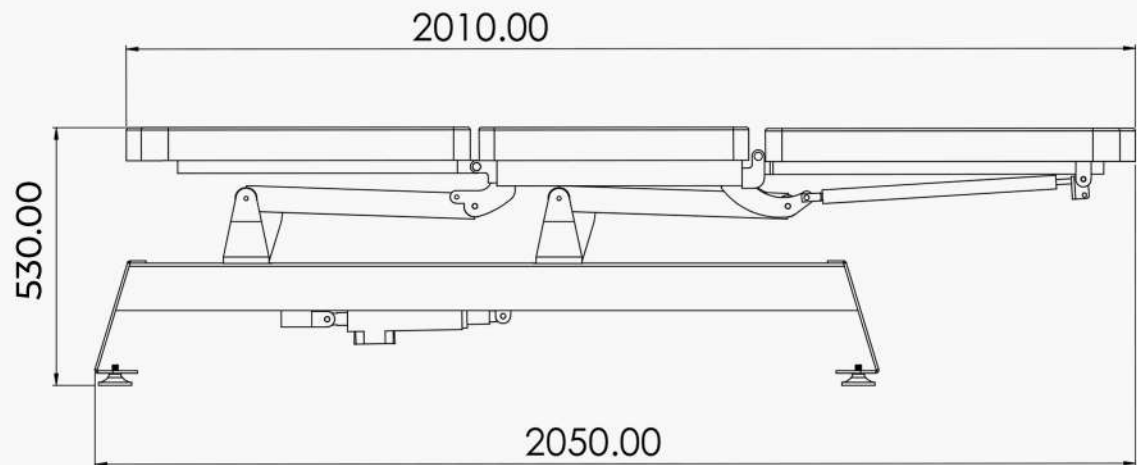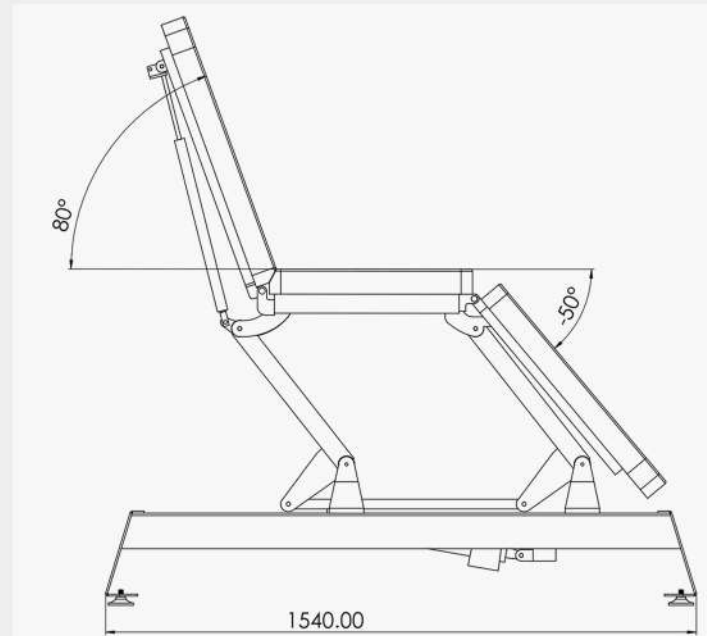


All dimensions are in mm



## MINIMUM HEIGHT

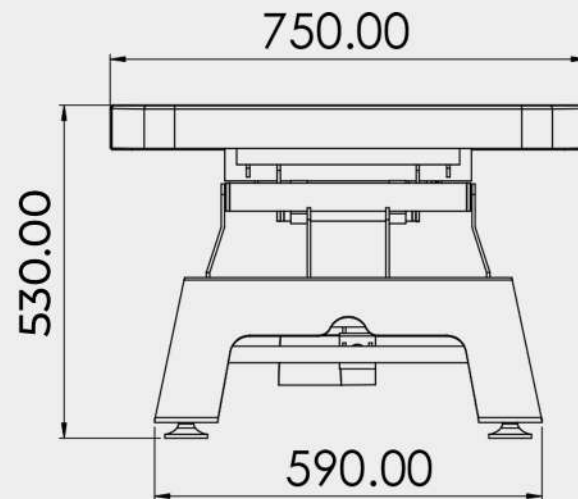

## MAXIMUM HEIGHT

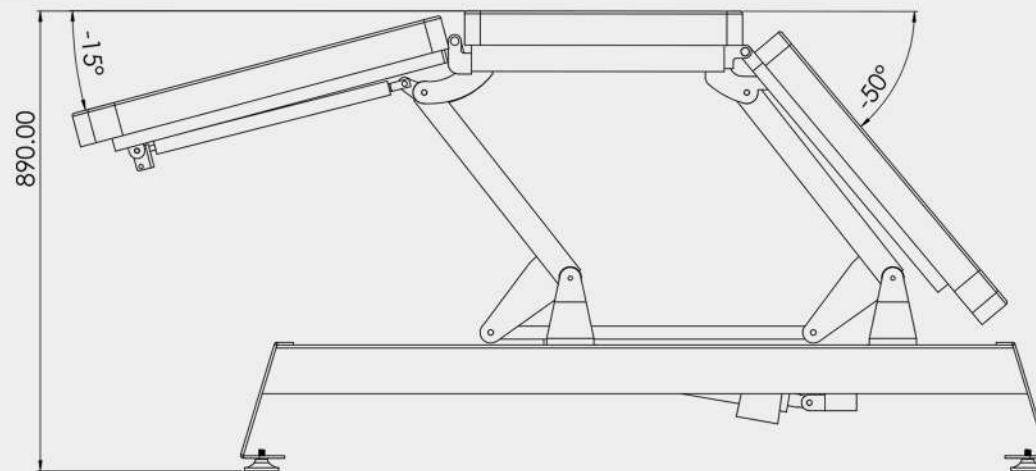

All dimensions are in mm

# MRIDUL

Derma Chair



OPTIONAL



In-Built Heating System

Professional 4 Motors Derma Chair with electronic height, backrest, leg rest and tilt function. This electronic dermatology chair comes with detachable headrest, detachable legrest, armrests, a nose slot for breathing and head pillow. The handheld remote can independently control height/backrest/ tilt/ leg rest position. The Lock function freezes the remote functions and the RESET button brings the Derma Chair back to its chair position.

# MRIDUL

01



02



03



04



**TOP:**

Top Design

Foam

Upholstery

Frame

Adjustment (Back Rest)

Adjustment (Height)

Adjustment (Leg Rest)

Five piece fully electronic Derma Chair with electronic height, backrest, leg rest and tilt function

75mm Mixed Layers for comfortable seating

PVC Coated Fabric

Mild Steel Construction Powder Coated  
Frame and ABS Molded Parts

Electronic

Electronic

Electronic

**BASE & FRAME:**

Base

Lifting Capacity

Powder Coated Metal Frame and  
ABS Components

Upto 360 Kgs

**UPHOLSTERY:**

Type

Available Colour Options

Weight


Abrasion

Flame/UV/Abrasion Resistant 1.2 mm  
Coated Fabric

25 colour options

660 gsm + 5%

EN ISO 5470=2, 25000 Cycles

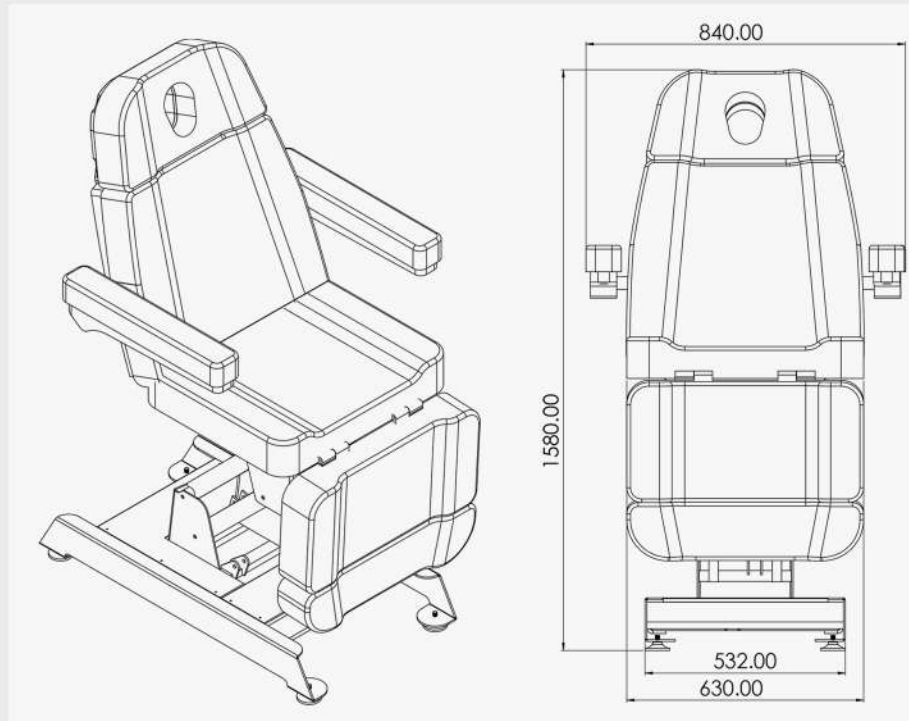**ACCESSORIES:**

Available Colour Options

25 colour options

**GUARANTEES:**

Structure

Electronic Motors &amp; Controls

1 year

3 years

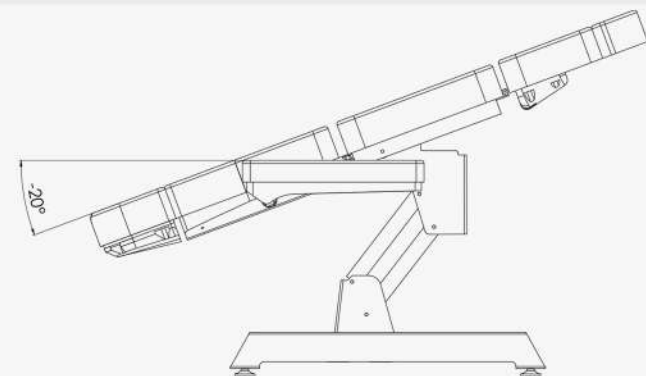
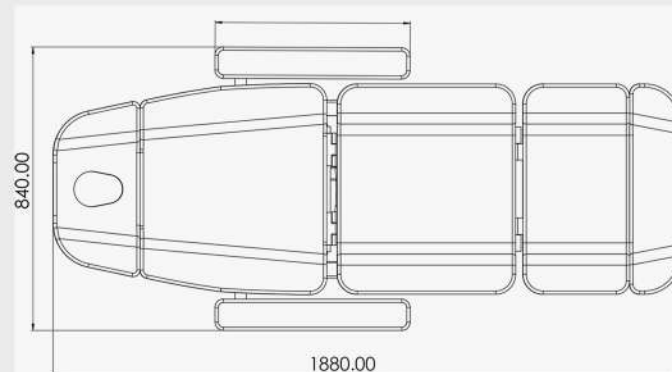
**WEIGHTS:**

Table Weight


Shipping Weight

90kgs


150kgs

**All dimensions are in mm**

## MINIMUM HEIGHT

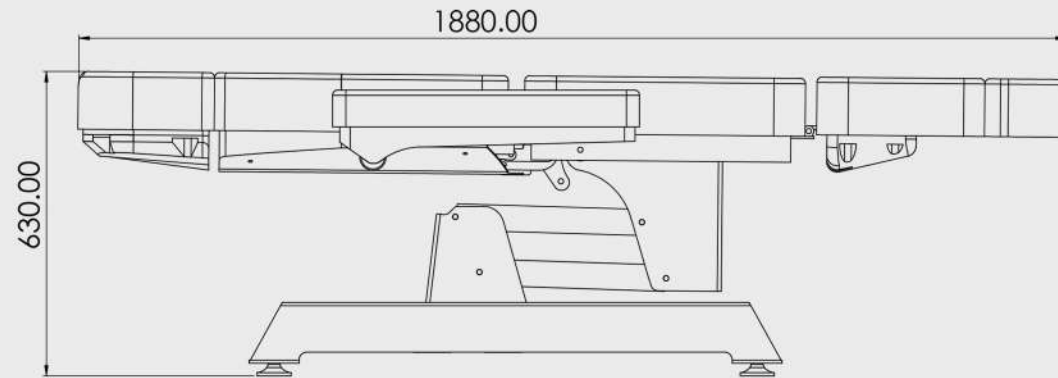

## MAXIMUM HEIGHT

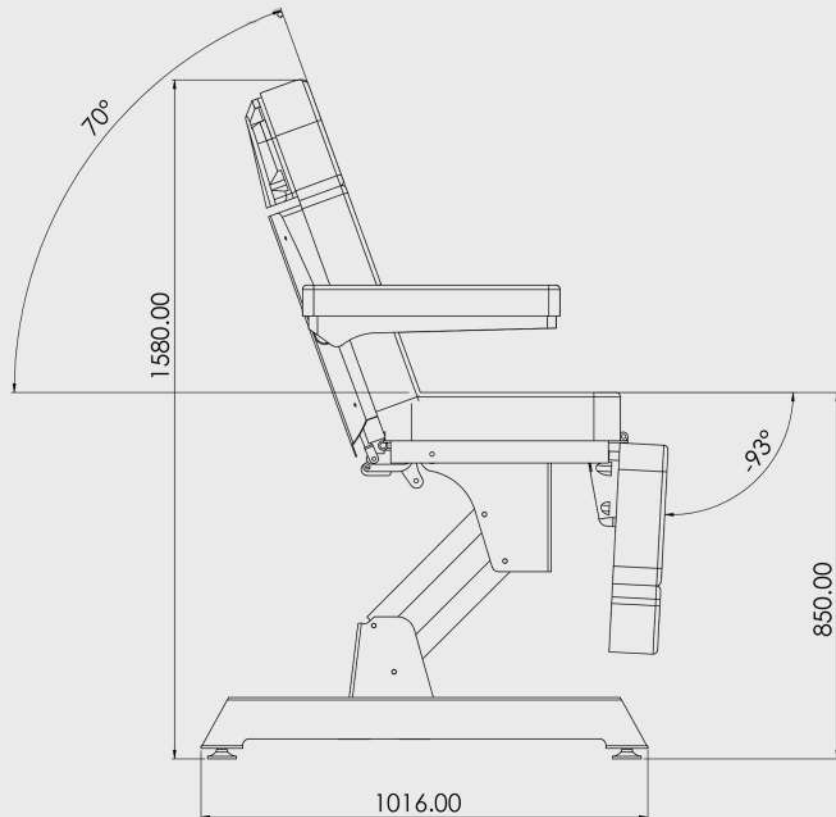

All dimensions are in mm



# GYNAE EXAMINATION CHAIR

Our Gynae Exam Chair, powered by three robust motors, offers seamless adjustment of height, backrest, and pelvis inclination for precise patient positioning. Crafted with patient comfort in mind, every detail, from the seat cushion to the footrests and headrest, ensures maximum comfort throughout the examination.



01



02



## TOP

Top Design	Four Piece articulating top
Top Cushion Measurements	1880 x 630 mm
Foam	75 mm mixed layer for comfortable seating
Upholstery	Nano Compound PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame & ABS Molded Parts

Adjustment (Back Rest)	Electronic
Adjustment (Leg Rest)	Detachable
Adjustment (Height)	Electronic
Adjustment (Tilt)	Electronic

## BASE & FRAME

Base	Powder Coated Metal Frame
Height Minimum	630 mm
Height Maximum	850 mm
Tilt Maximum	20 degrees
Lifting Capacity	Upto 360 kgs

## UPHOLSTERY




Type	Disinfectant Resistant / Chlorine Water Resistant / Abrasion / Flame / UV Resistant 1.10mm coated fabric with 2% Nano high-tech compound, 82% PV Comopound and 16% Knitted Polyester
Colour Options	25 colour options
Weight	820 gsm +/- 30 g/m2
Abrasion	ISO 12947-12 Kpa 2,00,000 Cycles

## ACCESSORIES

Detachable LED Light	Included with chair
Lithotomy Leg Rest	Included with chair


## GUARANTEES

Structure	5 year against any manufacturing defects
Upholstery	1 year against any manufacturing defects
Motors	3 years on Motors

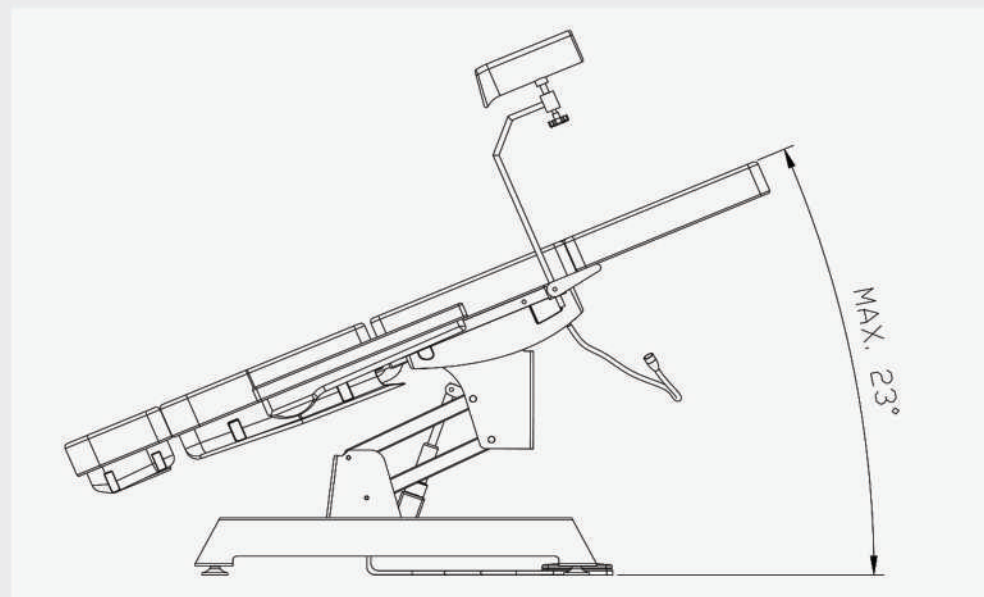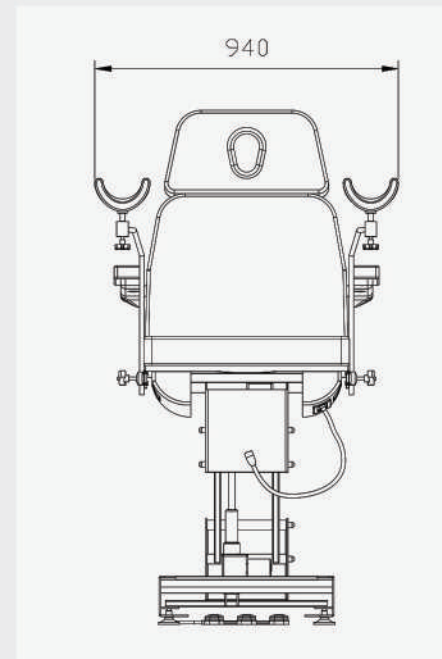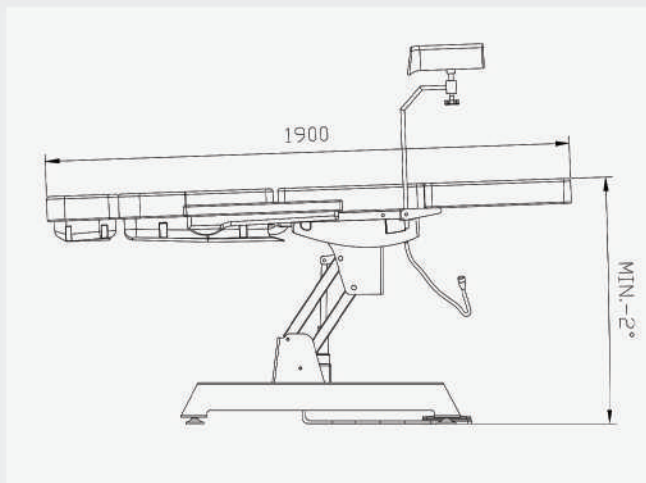


All dimensions are in mm

## MAXIMUM HEIGHT

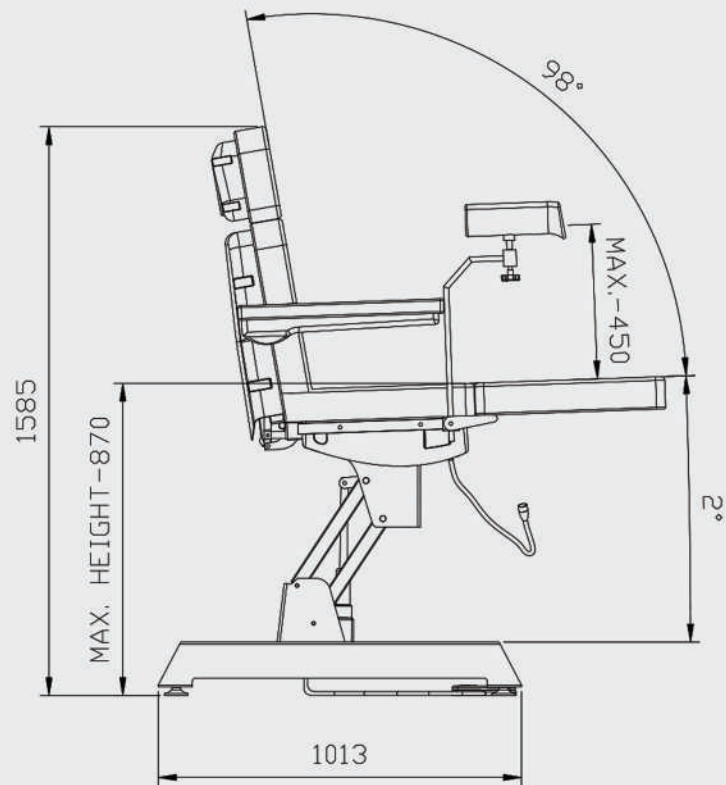

## MINIMUM HEIGHT

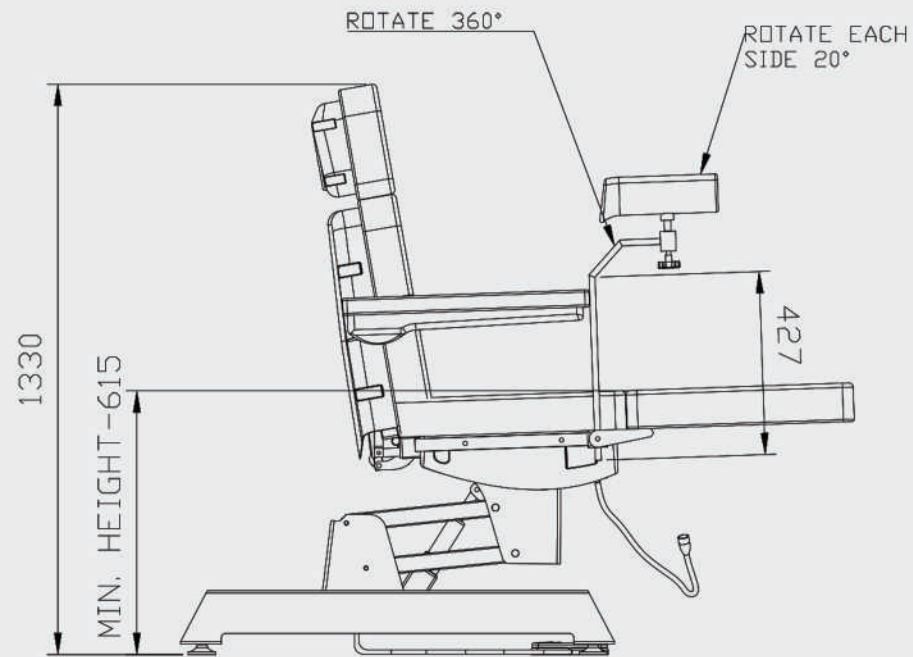

All dimensions are in mm



# MAHIT CHIROPRACTIC TABLE

Crafted with precision, this 7-section treatment table offers comprehensive adjustments for the head, back, legs, and mid-section of the plinth. Its electrically powered mechanism allows seamless elevation, granting clinicians the flexibility to manoeuvre into various positions for postural drainage and spinal manipulation with ease.





01



02





03



## TOP

Top Design	Seven Piece articulating top
Top Cushion Measurements	2020 x 760 mm
Foam	75 mm mixed layer for comfortable seating
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame & ABS Molded Parts
Adjustment (Back Rest)	Manual (Gas Spring)
Adjustment (Neck Rest)	Manual (Gas Spring)
Adjustment (Upper Back Rest)	Manual (Gas Spring)
Adjustment (Height)	Electronic (Linak Motor)

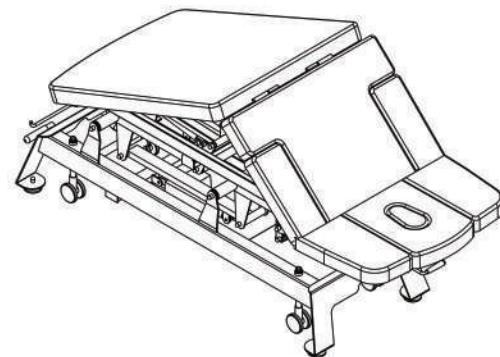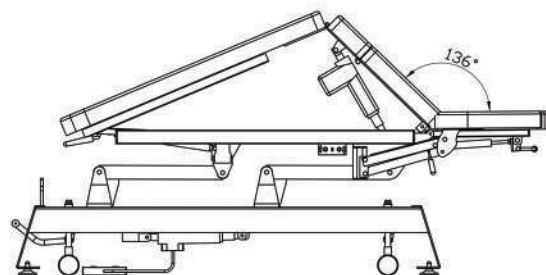

## BASE & FRAME

Base	Powder Coated Metal Frame
Height Minimum	530 mm
Height Maximum	1010 mm
Walk Free Switch	Enables you to operate the table from both sides
Lifting Capacity	Upto 220 kgs
Castor Wheels	For easy mobility in the treatment room

## UPHOLSTERY


Type	Abrasion / Flame / UV Resistant 1.2mm coated fabric
Colour Options	25 colour options
Weight	660 gsm + 5%
Abrasion	EN ISO 5470+2 25000 Cycles

## ACCESSORIES


Available Colour Options	25 colour options
--------------------------	-------------------

## GUARANTEES


Structure	5 year against any manufacturing defects
Upholstery	1 year against any manufacturing defects
Motors	4 years on Linak Motors

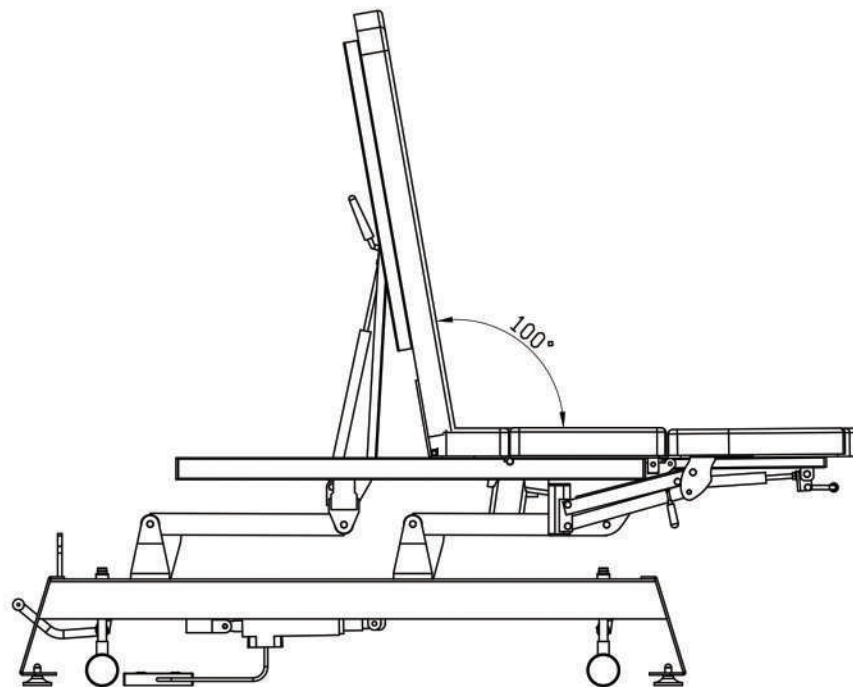

All dimensions are in mm

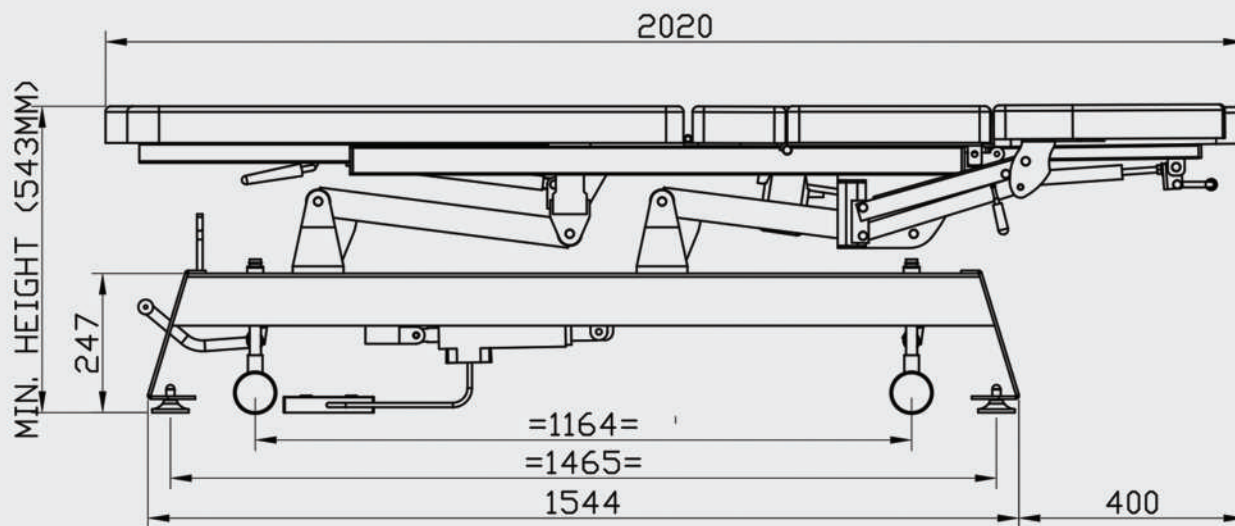

**MINIMUM HEIGHT**

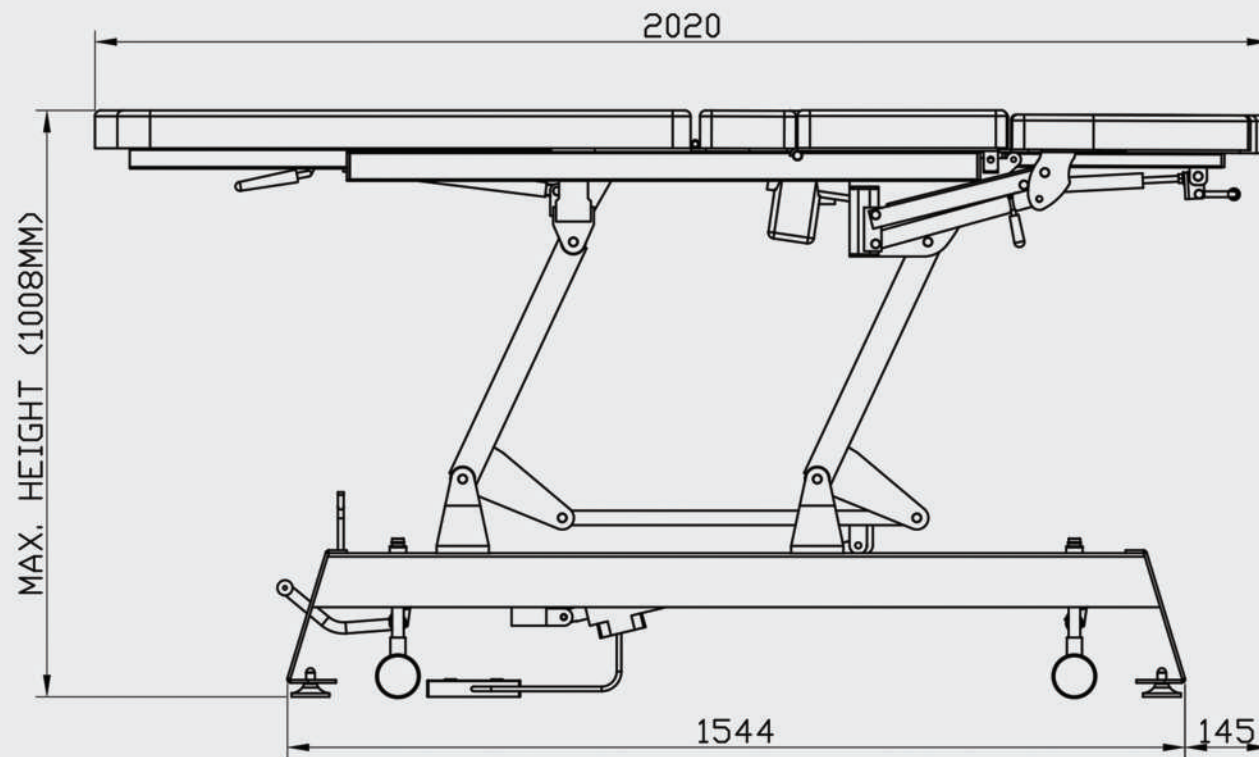

**MAXIMUM HEIGHT**

All dimensions are in mm



# DIALYSIS CHAIR

Our high-performance Dialysis Chair features three electric motors and lifting columns, seamlessly adjusting patients into their desired positions, including the emergency Trendelenburg stance, for safety and comfort. Designed for extended hemodialysis sessions, it boasts soft cushioning and disinfectant-resistant upholstery for impeccable hygiene.



01



02





## TOP

Top Design	Four Piece articulating top
Top Cushion Measurements	1890 x 560 mm (flat bed position) / 1540 x 560 (chair position)
Foam	75 mm mixed layer for comfortable seating
Upholstery	PVC Coated Fabric
Frame	Mild Steel Construction Powder Coated Frame & ABS Molded Parts
Adjustment (Back Rest)	Linak Denmark Quick Release Actuator
Adjustment (Leg Rest)	Linak Actuator
Adjustment (Height)	Linak Lifting Column
Adjustment (Armrest)	Detachable and adjustable armrests suitable for blood transfusion
Foot Rest	Adjustable/lockable foot pedestal to adjust according to patient's height

## BASE & FRAME

Base	Powder Coated Metal Frame
Frame Structure	Tubular frame with ABS Covers
Castor Wheels	Lockable Castor Wheels for easy mobility
Height Minimum	550 mm (flat bed position) / 1270 mm (chair position)
Height Maximum	850 mm (flat bed position) / 1560 mm (chair position)
Tilt Maximum	20 degrees
Lifting Capacity	Upto 360 kgs

## UPHOLSTERY


Type	Disinfectant Resistant / Chlorine Water Resistant / Abrasion / Flame / UV Resistant 1.10mm coated fabric with 2% Nano high-tech compound, 82% PV Comopound and 16% Knitted Polyester
Colour Options	25 colour options
Weight	820 gsm +/- 30 g/m2
Abrasion	ISO 12947-12 Kpa 2,00,000 Cycles

## ACCESSORIES

Head Cushion Support	Included with table
Adjustable phone & tablet holder	Included with table
Food Tray	Included with table
IV Stand Mount	Provision to attach an IV Stand


## GUARANTEES

Structure	5 year against any manufacturing defects
Upholstery	1 year against any manufacturing defects
Motors	4 years on Linak Motors

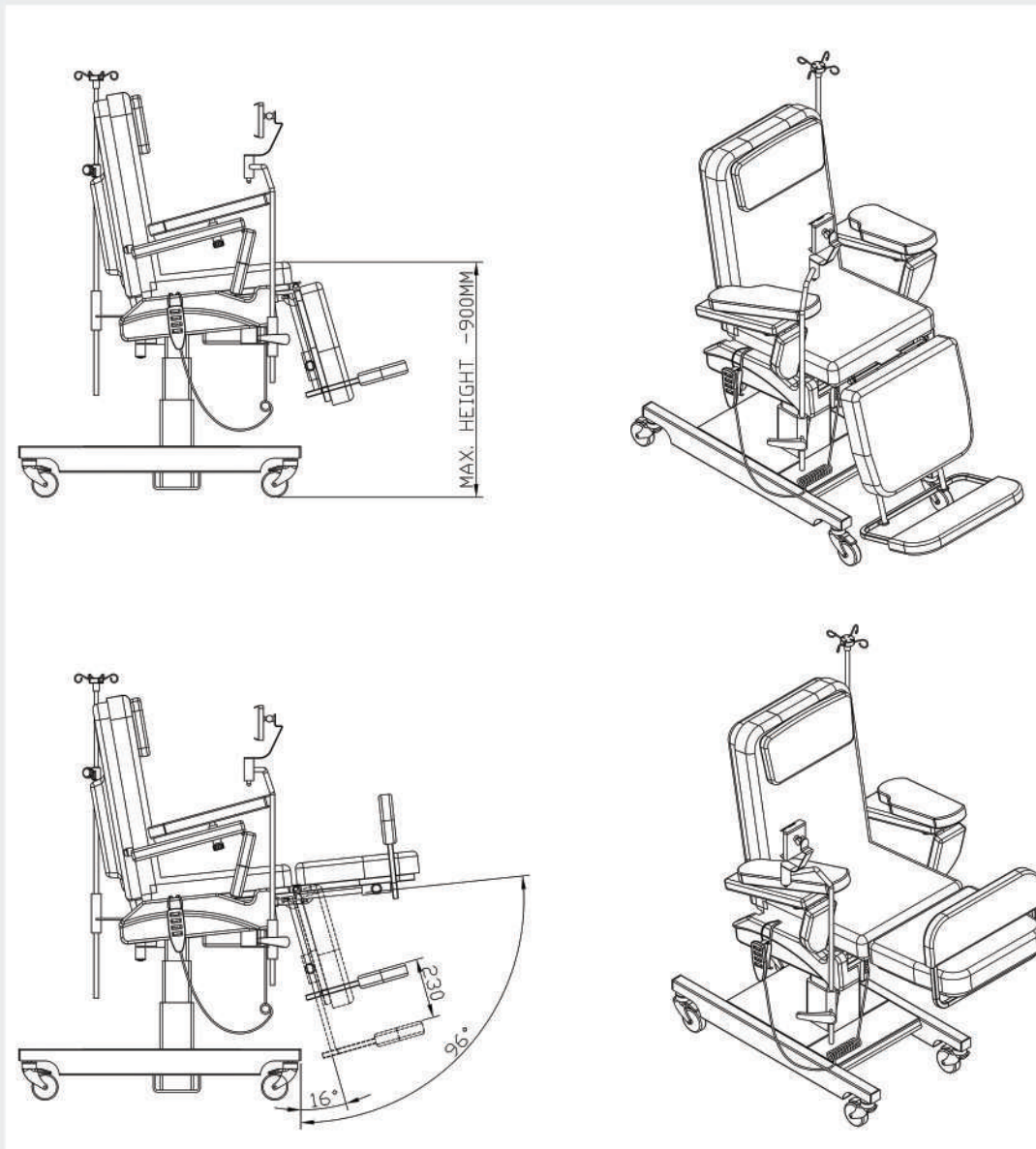


All dimensions are in mm

## MAXIMUM HEIGHT

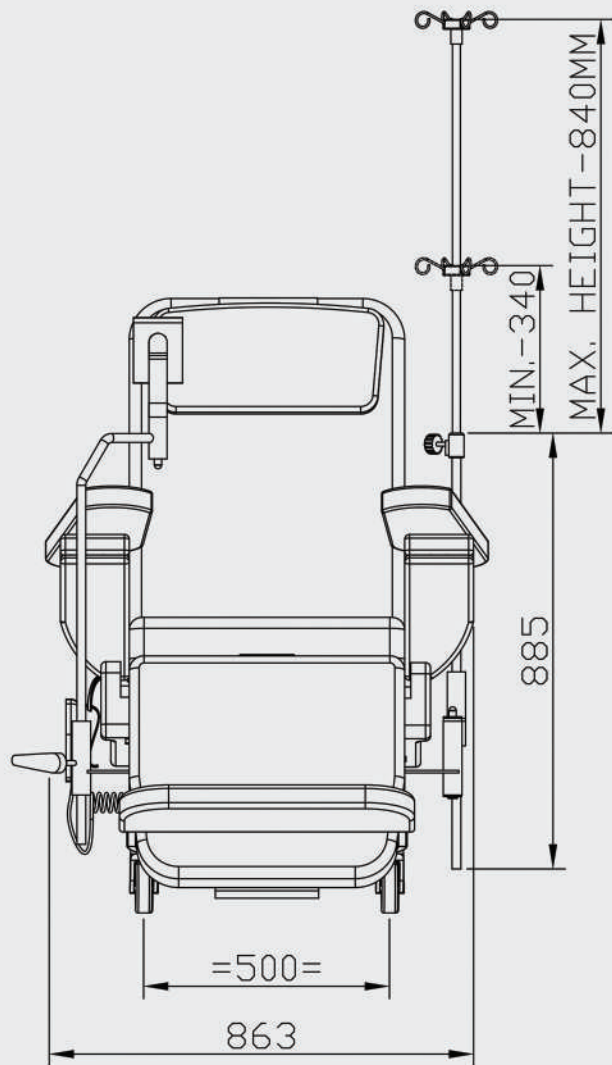

## MINIMUM HEIGHT

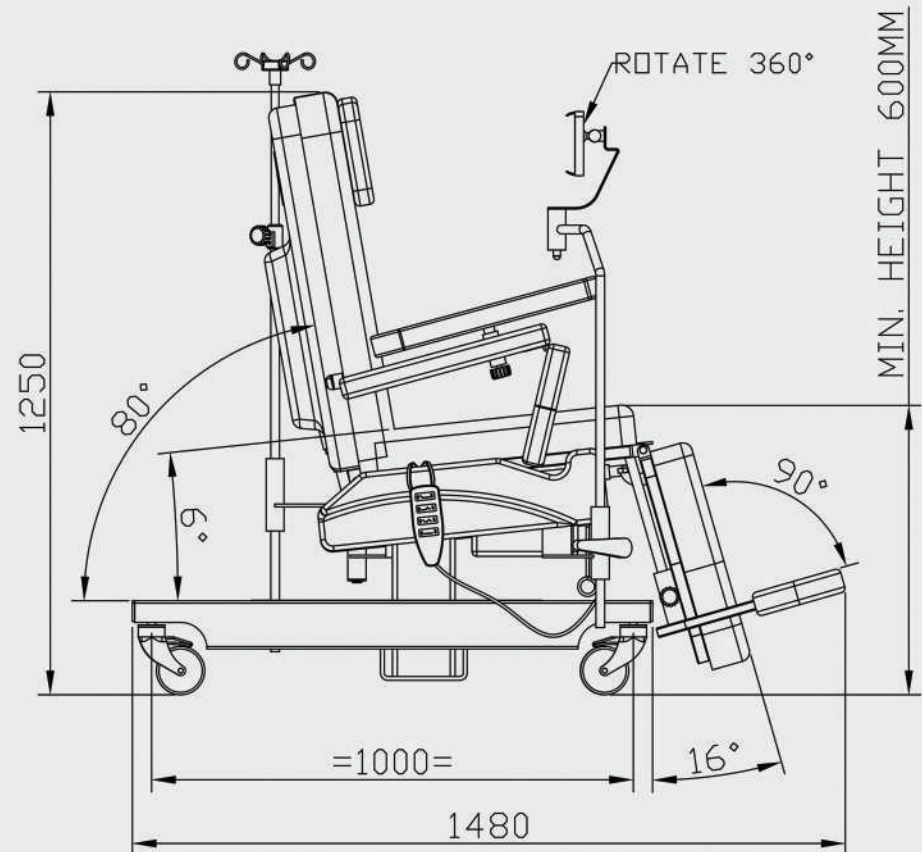

All dimensions are in mm